# Compat<sup>®</sup>

ENTERAL DELIVERY SYSTEM

## *Operating Instructions*

*For Model Number*

**9522041 / 199245** COMPAT SelectFlo<sup>®</sup>  
Enteral Delivery Pump  
and Memory

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## Table of Contents

COMPAT SelectFlo® Enteral Delivery Pump

### NOT FOR PARENTERAL USE

INTENDED EXCLUSIVELY FOR ENTERAL TUBE  
FEEDING FORMULAS.

### Section

Introduction	3
Controls	4
Operating Instructions	6
Alarms/Indicators	9
Troubleshooting	10
Specifications	11
Symbols on Pumps and Accessories	12
Self Test	13
Factory Service	13
Safety and Warnings	13
Limited Warranty	14
Cleaning	14
COMPAT® Enteral Delivery System	15

### For More Information

Contact your local Nestlé HealthCare Nutrition Sales  
Representative or call Customer Service from anywhere in the  
U.S. and Canada 1-877-338-4238.  
Monday–Friday 8:00 a.m.–4:30 p.m. Central Standard Time

Pump is not intended for neonatal use.

**Caution:** Federal law restricts this device to sale  
by or on the order of a licensed healthcare  
practitioner.

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
---

## Introduction

The COMPAT SelectFlo® Enteral Delivery Pump is an easy to use, portable pump for the administration of all enteral feeding solutions. Key features of this device include:

- **Latched Door**  
Reduces the likelihood of accidental dislodgement of sets and enables immediate detection of a free flow condition.
- **Automatic Priming**  
Frees up nurse to perform other bedside duties.
- **Ease of Use**  
Feeding parameters are easily displayed.
- **Versatility**  
Formula Rate adjustable in 1 mL/hr increments from 1–295.
- **Lightweight and Portable**  
The unit is 5.1 pounds, and can clamp to any IV pole.
- **Ease of Maintenance**  
The fire resistant plastic case and sealed control panel prevent fluid leakage and allow cleaning with any mild soap.
- **Easily Understood Alarms**  
Low Battery, Set Error and Flow Error for detecting empty container or occluded line.
- **2 ½ Minute Hold Alarm**
- **Low Occlusion Pressure**  
20 p.s.i. maximum occlusion allows for safe use with any enteral feeding tube.
- **Volume Delivered**  
For ease of tracking what has been administered when continuous feeding is desired, or over several feedings when intermittent feeding is necessary.
- **Memory Retention**  
of Infusion Rate and Volume Delivered even after unit is turned off. Settings are automatically recalled when the pump is turned back on.
- **Adjustable Alarm Volume**
- **Built-in Pole Clamp**
- **Luer Adapter Cap Holder**
- **Serviceability**  
Self-diagnostic mode allows for complete circuitry check when the pump is turned on.


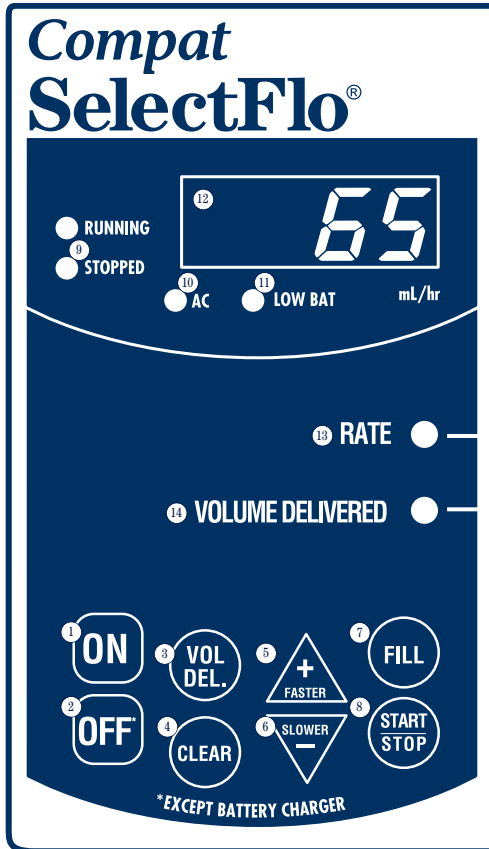
### Keypad Controls

1. **ON**, applies power to unit.
2. **OFF**, turns power off except to battery charger.
3. **VOL. DEL.**, press to display volume of formula delivered.
4. **CLEAR**, clears volume delivered to "0".
5. **▲** increases flow rate when pressed and held.
6. **▼** decreases flow rate when pressed and held.
7. **FILL**, automatically fills the administration set with formula.
8. **START/STOP**, alternately starts and stops the pump and silences the alarm.

### Function Indicator Displays

9. **RUNNING/STOPPED**, green RUNNING light indicates the pump is delivering formula. Red STOPPED light indicates the pump is stopped.
10. **AC**, indicates the unit is plugged into an AC outlet.
11. **LOW BAT**, indicates when the battery charge is below performance level.
12. **Display panel**, indicates selected parameter.
13. **RATE**, green RATE light indicates that the number in the display panel is the feeding rate.
14. **VOLUME DELIVERED**, green VOLUME DELIVERED light indicates that the number in the display is the amount of formula delivered.


## Controls (cont.)

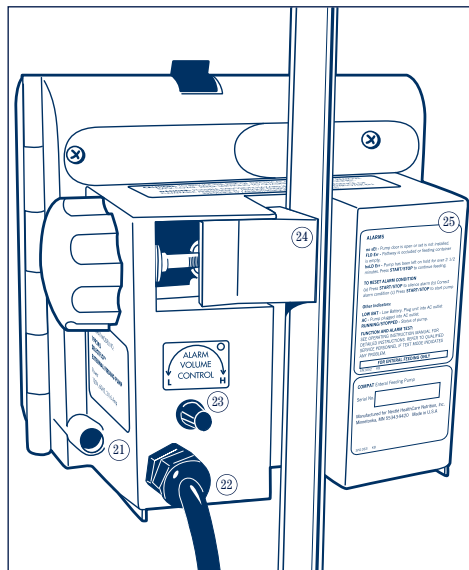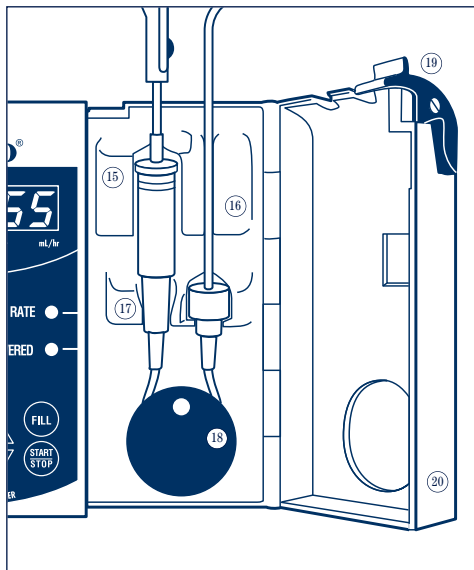15. Flow Sensors and Drip Chamber Guide
16. Tubing Guide
17. Drip Chamber Brackets
18. Rotary Peristaltic Pump Mechanism
19. Door Latch
20. Door
21. Luer Adapter Cap Holder
22. Power Cord
23. Alarm Volume Control
24. Pole Clamp
25. Directions (side panel also)


## Operating Instructions

6

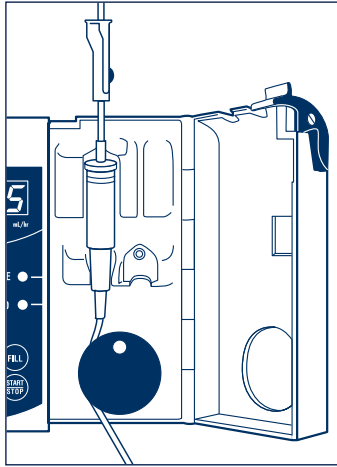
**DIAGRAM 1**

Open pump door by pushing, then pulling door latch, insert pump set drip chamber into drip chamber bracket.

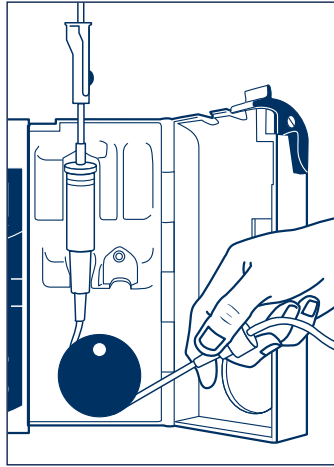
**DIAGRAM 2**

Assemble silicone tubing around rotor.

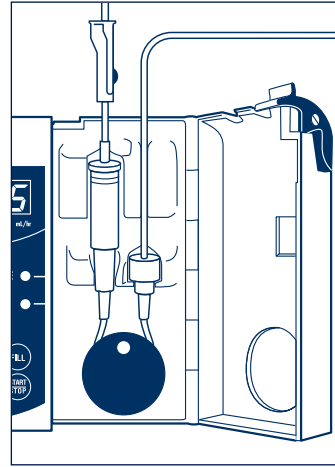**DIAGRAM 3**


Insert pump adapter into bridge and place tubing into tubing guide. Close door.

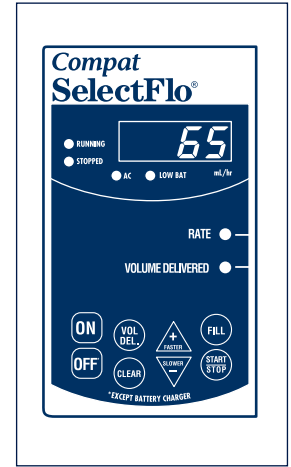**DIAGRAM 4**

Press ON to activate pump. Press ▲ (FASTER) or ▼ (SLOWER) until the desired setting is displayed. Clear Volume Delivered if necessary.

## Operating Instructions (cont.)

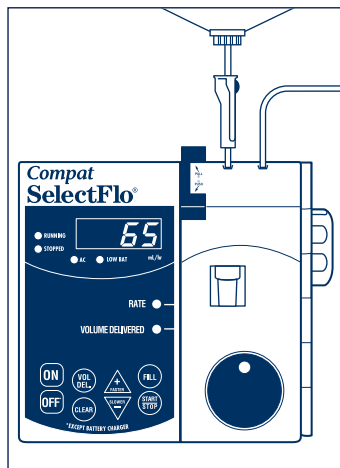**DIAGRAM 5**  
Completely assembled pump set.

Plug power cord into a 120 VAC, 60 Hz. three wire grounded outlet.

### Assemble Delivery Set

1. Open pump door by pushing/pulling door latch.
2. Prepare and hang the delivery set from an IV pole.
3. Insert delivery set drip chamber into drip chamber guide (roller clamp should be closed) (Diagram 1).
4. Assemble silicone tubing around the rotor (Diagram 2).
5. Insert pump adapter into bridge and place tubing into tubing guide (Diagram 3).
6. Close pump door.

### Check Pump Settings

1. Press ON—alarm will sound briefly and display '8888'. The pump has memory, and will scroll through the parameters and display the previous settings. It will return to flashing the Infusion Rate (Diagram 4).

2. Check the Infusion Rate. To change, press ▲ (FASTER) or ▼ (SLOWER) until desired rate is displayed.
3. Check Volume Delivered. Press VOL DEL. until Volume Delivered is indicated on the display panel. If desired, this can be reset to zero by pressing CLEAR.

### Auto Filling

1. Open the roller clamp on the SelectFlo delivery set. Remove the luer adapter cap from the delivery set, and place the cap in the luer adapter cap holder available on the pump.
2. Press the FILL key for automatic priming of the delivery set. The display will show the word FILL, and the rotor will begin to turn to indicate the unit is priming. Fill the entire line. The prime cycle takes approximately 2 minutes. The pump will audibly alarm and flash the infusion rate when filling is complete. (Note: Automatic filling may result in a small amount of formula being expelled beyond the luer adapter on the delivery set.)

### Gravity Fill

1. Follow steps 1–3 listed under **Assemble Delivery Set**.
2. Next, remove adapter cap from the delivery set, and place the cap in the adapter holder available on the pump. Open the roller clamp and fill the entire line with formula. Close the roller clamp.
3. Proceed to steps 4–6 under **Assemble Delivery Set**.
4. Complete steps 1–3 under **Check Pump Settings**.

### Begin Feeding

1. Connect the delivery set to the patient's feeding tube.
2. Press START/STOP to start the pump.
3. OFF will turn off the pump at any time (Infusion Rate and Volume Delivered are retained).

---

## **Operating Instructions (cont.)**

### **How to Modify Pump Settings**

#### **To change the Infusion Rate**

1. Press START/STOP to stop the pump.
2. Enter a new infusion rate by pressing ▲ (increase) or ▼ (decrease) until the desired rate is displayed.

#### **To clear Volume Delivered**

Volume Delivered is primarily intended to monitor formula delivered when intermittent feedings are required.

1. Press START/STOP to stop the pump if running.
2. Press VOL DEL until Volume Delivered light is lit.
3. Press CLEAR.

---

## Alarms/Indicators

### Note

**Pressing the START/STOP key after the alarm condition is identified will silence the alarm and return pump to hold mode.**

### AC

Indicates the pump is operating on electrical power. The pump should be left plugged into an AC outlet whenever practical. No audible alarm will sound.

### Running/Stopped

Indicates the status of the pump.

### Low Battery (LOW BAT)

This is activated when the battery is run down. After 15 minutes, this light will flash, an audible alarm will sound and the pump will stop. Whenever the LOW BAT light is on, the pump should be plugged into an AC outlet, and normal feeding may continue. To avoid degrading battery life it is important not to run down the battery. The battery recharges anytime the pump is connected to an AC outlet, regardless of whether the pump is running.

### Flow Error (FLO Err)

Under these conditions the pump stops, an alarm sounds and is displayed as 'FLO Err'. One of the following conditions may apply:

#### Possible Condition

1. Feeding container is empty
2. An occlusion has caused delivery to stop
3. Drip chamber is overfilled
4. Drip chamber walls are severely coated with formula
5. Sensors in drip chamber bracket need cleaning
6. Drip chamber is improperly placed

#### Corrective Action

1. Refill container or start new administration set
2. Identify point of occlusion
3. Insert new administration set
4. Manipulate chamber to clear drops in sensor path
5. Clean sensor with isopropyl alcohol
6. Reposition drip chamber

---

### No Set (no SET)

If the door is not closed and latched securely, no set is loaded on the pump, or a proper COMPAT SelectFlo® set is not used, the alarm message 'no SET' will be displayed. Correct by loading an appropriate COMPAT SelectFlo delivery set, close and latch the door, and press START/STOP.

### Hold error (hoLd Err)

If pump is left on hold and START/STOP is not pressed within 2 1/2 minutes, this alarm will be activated and 'hoLd Err' will be displayed. Press START/STOP button twice to extend the hold for another 2 1/2 minutes; or press once to restart the feeding pump.


## Troubleshooting

10

COMPAT® Enteral Delivery Pumps are reliable electromechanical devices. As with any mechanical device, problems may occur. The following tips will help you correct the minor in-service problems that may arise.

TROUBLESHOOTING GUIDE		
SITUATION	POSSIBLE CAUSE	SOLUTION
<ul style="list-style-type: none"> <li>• Low battery alarm activated</li> <li>• (LOW BAT) illuminated</li> </ul>	• Battery charge is below performance level	• Plug cord into AC outlet
	• Feeding container is empty	• Replace or refill feeding container
<ul style="list-style-type: none"> <li>• (FLO Err) alarm activated</li> </ul>	• Occlusion is restricting flow	<ul style="list-style-type: none"> <li>• Identify point of occlusion and free up               <ul style="list-style-type: none"> <li>- Pump Set                   <ul style="list-style-type: none"> <li>• kinked tubing</li> <li>• blocked tubing</li> <li>• closed roller clamp</li> </ul> </li> <li>- Feeding tube</li> <li>- Feeding container</li> </ul> </li> </ul>
	• Drip chamber not properly placed in pump	• Reposition drip chamber
	• Drip chamber walls coated with formula drops or condensation in sensor path	• Manipulate chamber to clear drops or condensation in sensor path
	• Flow sensors housing need cleaning	• Clean with cotton swab dampened with isopropyl alcohol
	• (hold Err) Pump is left on hold	• Press START/STOP button
<ul style="list-style-type: none"> <li>• (No Set) Alarm activated</li> </ul>	• Door not closed and latched	• Examine tubing for proper placement, close door
	• Wrong set is used	• Load an appropriate set
	• Set is not properly in place	• Assemble tubing properly around pump rotor. Close and latch door

Specifications

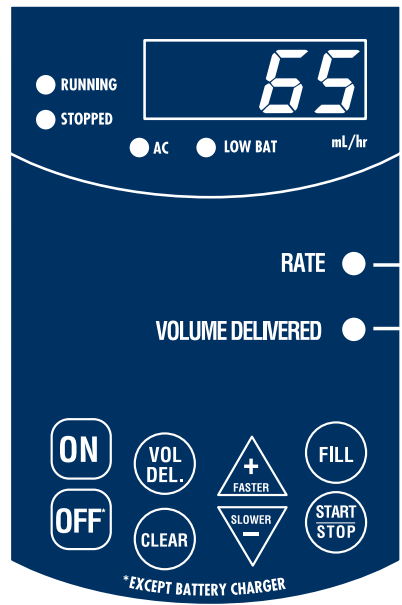Type of Device .....	Volumetric enteral feeding pump
Mechanism .....	Rotary Peristaltic
Infusion Rate .....	1–295 mL/hr in 1 mL increments
Accuracy .....	± 10% of selected flow rate
Occlusion Pressure .....	Does not exceed 20 p.s.i.
Battery .....	Sealed lead acid Rechargeable Automatically recharges when connected to an AC outlet Recharge time is 12 hours from complete discharge to full recharge Operating time is 10 hours at 100 mL/hour (fully charged battery) (Charge unit fully every 6 months)
Dimension .....	6.4" High x 6.4" Wide x 5.0" Deep
Weight .....	5.1 pounds
Case Material .....	Fire and impact resistant plastic
Power Requirement .....	120V~, 187 mA, 60 Hz
Operating Temperature .....	10°–40°C (50°–104°F) 75% R.H. non-condensing
Storage Temperature .....	0°–50°C (32°–122°F) 95% R.H. non-condensing
Type of Protection Against Electrical Shock .....	Class 1, Internally Powered Equipment
Degree of Protection Against Electrical Shock .....	Type BF
Degree of Protection Against Ingress of Fluids .....	Ordinary Equipment – IPX0

---

## Symbols on Pump and Accessories

12

The following symbols may be found on the pump:

Refer to Accompanying Documents

Alternating Current

Type BF Protection (degree of protection against electrical shock—there is no conductive connection to the patient)

Medical Electrical Equipment

- 1) Classified with respect to electrical shock, fire and mechanical hazards in accordance with UL60601-1
- 2) Classified with respect to electrical shock, fire, mechanical and other specified hazards in accordance with CAN/CSA C22.2 No. 601.1

---

## Self Test

The major functions of the pump may be tested by the following steps:

1. Turn the pump on by pressing the ON key.
2. 8888 will be displayed immediately. While this is displayed, press and simultaneously hold the ▲ and ▼ keys for 4 seconds.
3. After 4 seconds, an audible alarm will sound, and the screen will display 'teSt' momentarily, and then all eights (8888).
4. The pump rotor should begin to turn. To test the visual alarms, press either CLEAR, ▲, ▼, FILL, VOL DEL, or START/STOP. When any of these keys are pressed or the sensor is momentarily blocked, all visual alarm indicators should light, an audible alarm will sound, and the rotor will halt momentarily.
5. Any failure indicates a service check should be performed. See the Factory Service section for instructions on where to send the pump.
6. Turn the pump off by pressing the OFF key.

---

## Factory Service

In the event a pump requires servicing, contact your Nestlé HealthCare Nutrition Sales Representative or call Customer

Service for instructions from anywhere in the U.S. and Canada at 1-877-338-4238. All service work should be performed by authorized

personnel only. Unauthorized repairs void warranty (see section on limited warranty).

---

## Safety and Warnings

### **NOT FOR PARENTERAL USE INTENDED EXCLUSIVELY FOR ENTERAL TUBE FEEDING FORMULAS.**

Do not use this device near flammable anesthetics, it is a

possible explosion hazard.

This device is designed to minimize the effects of uncontrolled electromagnetic interference and other types of interference from external sources. Avoid use of other equipment that may cause erratic operation or degradation in the perfor-

mance of this device.

If the integrity of the external protective conductor in the installation is in doubt, equipment shall be operated from its internal electrical power source.

---

## Limited Warranty

14

Nestlé HealthCare Nutrition warrants the COMPAT® Enteral Feeding Pump against defects in material and workmanship under normal use and service for a period of one (1) year from date of delivery. This warranty is valid only to the original purchaser and does not extend to any product or part thereof which

has been subjected to accident, abuse, alteration, misuse, or which has not been operated and maintained in accordance with prescribed instructions. This warranty shall not apply if the product has been repaired by anyone other than an authorized service representative of Nestlé HealthCare Nutrition.

Nestlé HealthCare Nutrition reserves the right to repair or replace (at its option) any pump which fails to meet the foregoing warranty.

This warranty is in lieu of all other warranties expressed or implied including, without limitation, warranties of merchantability and of fitness for a particular purpose.

---

## Cleaning

### Caution

**Always disconnect the pump from the AC outlet before cleaning to avoid electric shock hazard.**

Using a cloth or sponge, clean the pump housing and rotor regularly with warm soapy water. (Prolonged exposure to alcohol, household detergents or strong cleansers can result in damage to pump housing.)

Clean drip sensors with cotton swabs and isopropyl alcohol.

### Disinfecting

If it is necessary to disinfect the entire pump, the following is recommended:

#### General Disinfectants

Cidex®  
pHisohex®  
Hibiclens®

Isopropyl alcohol

#### AIDS and Hepatitis Disinfectants

10% concentration of 5.25% sodium hypochlorite (household bleach)

#### Tuberculosis Disinfectants

70% concentration of  
isopropyl alcohol

After application of these disinfecting agents, allow to air dry. It should be noted that these recommendations should not be substituted for official procedures that may differ among institutions. To determine the proper procedure in a particular institution, contact the Infection Control Department.

**Do not autoclave. Do not immerse pump in cleaning solutions. Delivery sets are for single use only. Do not clean and reuse.**

Cidex®, pHisohex® and Hibiclens® are trademarks of their respective companies

---

## COMPAT® Enteral Delivery System

### COMPAT® Pumps

9517455 / 199235	Compat® Enteral Delivery Pump with Dose Limit
9522041 / 199245	Compat SelectFlo® Portable Pump featuring Latched Door
9521437 / 199255	Compat DualFlo® Enteral Delivery Pump with Hydration Feature

### COMPAT® Delivery Sets

#### Gravity Sets

9516074 / 199216	1000 mL Twist Cap Vinyl Bag with Gravity Set
------------------	--

#### For use with the 9517455 / 199235 Compat® pump and 9522041 / 199245 Compat SelectFlo® pump

9521433 / 199507	Pump Set with in line “Y” Adapter and Piercing Spike
9519810 / 199547	Pump Set with in line “Y” Adapter and 1000 mL Vinyl Bag
9519805 / 199312	500 mL Baggle Semi-Rigid Container with Pump Set
9522610 / 199707	SpikeRight™ Piercing Spike with in line “Y” Adapter and preattached Pump Set

#### For use with the 9521437 / 199255 Compat DualFlo® pump

9521429 / 199407	DualFlo Piercing Spike and 1000 mL Vinyl Water Bag
9521430 / 199447	DualFlo 1000 mL Formula Vinyl Bag and 1000 mL Vinyl Water Bag
9521434 / 199607	DualFlo with SpikeRight™ Piercing Spike and 1000mL Vinyl Water Bag

### For More Information

Call your Nestlé HealthCare Nutrition Sales Representative or Customer Service at 1-877-338-4238.

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Except as noted, trademarks are owned by
Nestlé HealthCare Nutrition, Inc.
Société des Produits Nestlé, S.A., Vevey, Switzerland.